



autotilt **sigma**



Made in
Great Britain

The specialist chair for critical
and bariatric care

The specialist chair for critical

The new range of battery powered Sigmas, which are designed to meet patient needs in Critical Care units, are ideal for enabling staff to move a patient out of bed at the earliest opportunity to aid recovery. From the operator's perspective, the powered controls make for easier movement and optimum patient positioning. In addition the Quick Release function gives staff the peace of mind that they can bring a patient into a horizontal position in under 3 seconds.

There is a choice of models for patients small or large, slight or obese.



MAIN FEATURES

- Motorised chair and trolley modes.
- Trendelenburg tilting $\pm 10^\circ$.
- Anti-bacterial upholstery and paintwork.
- Easily removable snap-lock safety belt.
- Castors can be set to steer front or rear, free movement, or locked in lateral pairs.



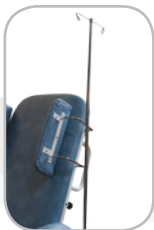
al care

autotilt **sigma****TILT IN SPACE**

The battery motors operate the variable height, chair to trolley mode and Trendelenburg/Reverse Trendelenburg positioning; $\pm 10^\circ$. The leg rests are independently adjustable for angulation and length.

**QUICK RELEASE**

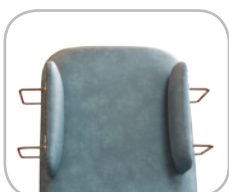
Manual quick release facility for horizontal positioning in seconds.

**Accessories****DRIP STAND**

- Double hook.
- 2 section adjustable height.
- Can be fitted to either side of the Sigma.

**OXYGEN BOTTLE HOLDER**

- Fits white CD2 O₂ cylinder.
- Can be fitted vertically to either side of the Sigma.

**LATERAL SUPPORTS (PAIR)**

- Adjustable inboard reach.
- Matching anti-bacterial upholstery to the Sigma.
- 2 positions: Upper and lower or backrest.



Specifications



autotilt **sigma 135**

Maximum Patient Weight Capacity	135 kg / 21 stone
Chair Width / Depth (Inside)	54 cm / 44 cm
Chair Mode Height Range	65 - 86 cm
Trolley Mode Height Range	84 - 106 cm



autotilt **sigma 250**

Maximum Patient Weight Capacity	250 kg / 39 stone
Chair Width / Depth (Inside)	74 cm / 44 cm
Chair Mode Height Range	64 - 84 cm
Trolley Mode Height Range	82 - 103 cm



autotilt **sigma 300**

Maximum Patient Weight Capacity	300 kg / 47 stone
Chair Width / Depth (Inside)	74 cm / 54 cm
Chair Mode Height Range	66 - 86 cm
Trolley Mode Height Range	84 - 106 cm