



Diagnostic Imaging Tables

www.knightimaging.co.uk



**Knight
Imaging**

EXPERIENCE & EXPERTISE

Imaging and C-Arm Tables

Multi-purpose Examination Table	3
C-ARM Tables	3

Specialist Tables

Shift CT Trolley	4
Motorised X-ray Trolley	4
Motorised Side-Access C-Arm Trolley	5

CT Tables

ACITOS (Accident Total Scan) CT Trolley SE	6
ACITOS (Accident Total Scan) CT Trolley EE	6

MRI Tables

MRI Trolley EE	7
MRI Trolley SS	7
MRI Trolley SSM	7


Multi-purpose Examination Table

Ideal for use in endoscopy, colonoscopy, radiology and pain clinics.

- Battery powered – no trailing wires.
- Battery driven motorised height adjustment from 60 cm to 85 cm.
- Battery driven motorised $\pm 10^\circ$ Trendelenburg.
- Floating X-ray transparent tabletop with electromagnet locks.
- X-ray transparent area: 64 cm x 197 cm.
- Flexible control panel on side rail.
- Extra battery with wall charger.
- Patient weight up to 250 kg.

Optional Fittings

- Carbon fibre tabletop.
- Side rails.
- Drip stand.
- Mattress with positioning handles available.
- Available in a choice of upholstery colours.

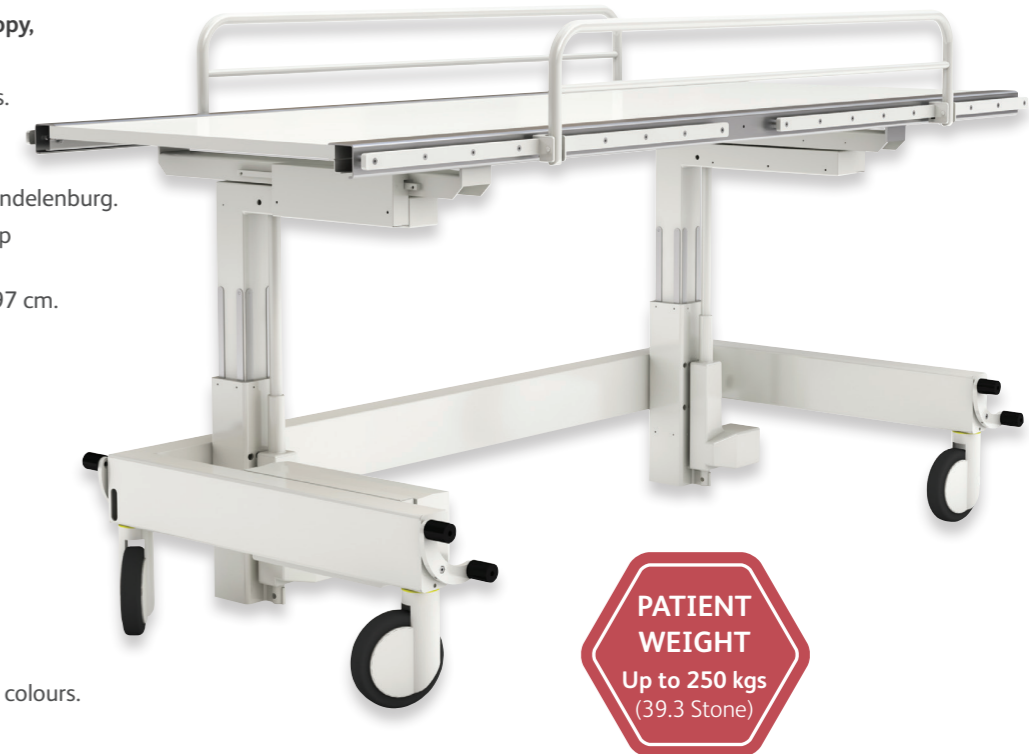

C-Arm Tables

A bespoke range of tables, available with up to 5 ranges of movement.

- Battery powered – no trailing wires.
- Battery driven motorised programmable working height from 70 to 100 cm.
- Carbon fibre X-ray transparent table with head profile: 56 x 160 cm.
- Wheels with electrical lock.
- Extra battery with wall charger.
- Patient weight up to 250 kg.

Optional Fittings

- Battery driven motorised 20 cm transversal and 40 cm longitudinal movement.
- Battery driven motorised $\pm 15^\circ$ Trendelenburg and $\pm 15^\circ$ transversal tilting.
- Joystick panel for precise control of tabletop movements.




Shift CT Trolley

Designed for maximum flexibility allowing all round imaging access. Operatives can switch effortlessly between imaging techniques without moving the patient.

The detachable carbon fibre tabletop reduces the moving and handling risk as it allows the patient to be scanned comfortably whilst remaining on the table.

- Fully detachable tabletop.
- Battery powered – no trailing wires.
- Detachable carbon fibre table with manual table locks: 56 x 205 cm.
- Battery driven motorised transversal movement of table top: 50 cm.
- Battery driven motorised transversal movement of columns: 50 cm.
- Battery driven motorised up/down movement of table top: 10 cm.
- Extra battery with wall charger.
- Patient weight up to 250 kg.

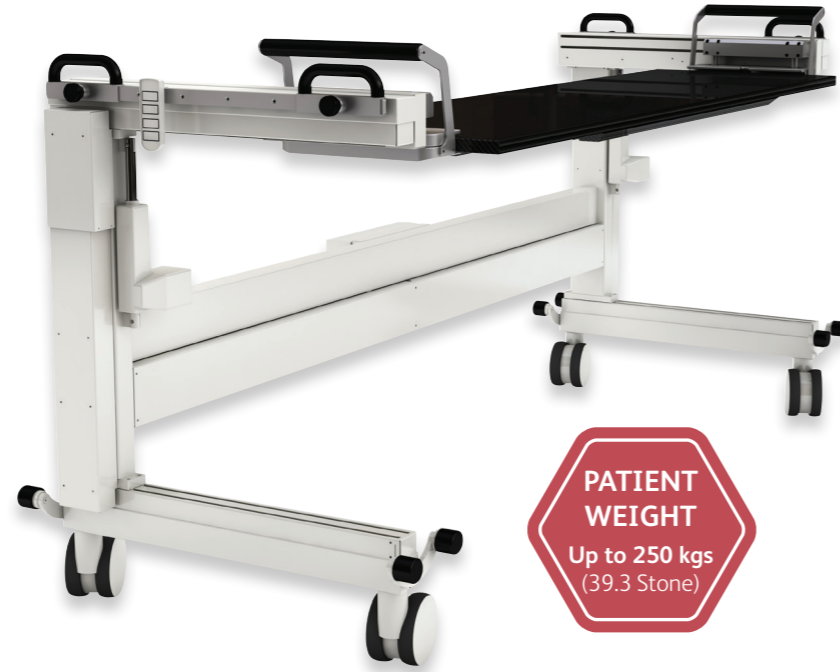

PATIENT WEIGHT
Up to 250 kgs
(39.3 Stone)

Motorised Side-Access C-Arm Trolley

Allows perfect positioning of a patient for X-raying with a C-Arm unit.

- Battery powered – no trailing wires.
- Tabletop with X-ray transparent area: 56 cm x 200 cm.
- Battery driven motorised transversal movement: 20 cm.
- Battery driven motorised longitudinal movement: 30 cm.
- Battery driven motorised working heights from 60 cm to 90 cm.
- Battery driven motorised $\pm 8^\circ$ Trendelenburg.
- 15 cm dual castors with electrical brakes and 2 control panels.
- Extra battery with wall charger.
- Patient weight up to 250 kg.

Optional Fittings

- Carbon fibre tabletop.
- Joystick panel for precise control of tabletop movements.



PATIENT WEIGHT
Up to 250 kgs
(39.3 Stone)

Motorised X-ray Trolley

Allows X-raying of a patient with a C-Arm unit or detector unit.

- Battery powered – no trailing wires.
- Fixed tabletop with X-ray transparent area: 56 cm x 200 cm.
- Battery driven motorised working height from 65 to 95 cm.
- 3 programmable working heights.
- Battery driven motorised $\pm 12^\circ$ Trendelenburg.
- Extra battery with wall charger.
- Large wheels with total and directional lock.
- Patient weight up to 250 kg.

Optional Fittings

- Eltronax tabletop with X-ray transparent area: 188 x 42 cm.



PATIENT WEIGHT
Up to 250 kgs
(39.3 Stone)

ACTIOS (Accident Total Scan) CT Trolley SE

For transportation and delivery of a patient directly to the CT scanner. The detachable carbon fibre tabletop reduces the moving and handling risk as it allows the patient to be scanned comfortably whilst remaining on the table.

- Battery powered – fully portable and no trip risk.
- Detachable carbon fibre table with manual table locks: 205 x 56 cm.
- Motorised height from 65 to 95 cm.
- 3 programmable heights.
- Motorised Trendelenburg.
- 15 cm wheels with central and total lock.
- Additional back-up battery with wall charger.
- Weight capacity up to 250 kg.



ACTIOS (Accident Total Scan) CT Trolley EE

For transportation and delivery of a patient from an accident and emergency environment directly to the CT scanner. The detachable tabletop allows the patient to be scanned whilst remaining in a stationary position.

- Battery powered – fully portable and no trip risk.
- Detachable carbon fibre table with manual table locks: 205 x 56 cm.
- Motorised height from 75 to 105 cm.
- 3 programmable heights.
- Outer width: 92 cm. Inner width (between arms): 70 cm.
- 15 cm wheels with central and total lock.
- Additional back-up battery with wall charger.
- Weight capacity up to 250 kg.
- Extra battery with wall charger.
- Patient weight up to 200 kg.



MRI Trolley EE

Ideal for delivering a patient to the MRI scanner. The detachable tabletop allows the patient to be scanned whilst remaining on the table.

- MRI safe – 3T.
- Detachable tabletop: 205 x 56 cm.
- Outer width: 92 cm. Inner width (between arms): 70 cm.
- Hydraulic pump adjustment of table height from 75 to 105 cm.
- Large wheels with central and total lock.
- Weight capacity up to 250 kg.



MRI Trolley SS

Delivers a patient to the MRI scanner to be scanned with a mattress.

- MRI safe – 3T.
- Tabletop with plastic side rail protection: 65 x 210 cm.
- Hydraulic pedal adjustment of table height from 65 to 95 cm.
- Large wheels with central and total lock.
- Extra battery with wall charger.
- Patient weight up to 200 kg.



MRI Trolley SSM

Delivers a patient to the MRI scanner to be scanned with a mattress.

- MRI safe – 3T.
- Tabletop (including mattress with 4 handles at each side): 56 x 200 cm.
- Hydraulic pedal adjustment of table height from 60 to 90 cm.
- Adjustable back rest and fold down side guard.
- Large wheels with central and total lock.
- Extra battery with wall charger.
- Patient weight up to 200 kg.



Service and Support

Knight Imaging offers service and support on our range of couches and chairs. We carry out one-off service calls, yearly preventative maintenance visits as well as offering support contracts.

Planned Preventative Maintenance

Our fully trained engineers service your couches and chairs ensuring that any potential issues are spotted and brought to your attention before they become a problem.

We typically carry out these types of services across entire departments ensuring the best possible lifetime of your couches and chairs. We will always arrange visits in advance to ensure a swift turnaround and to minimise disruption to your ward or department.

A preventative maintenance visit can be arranged whenever you think you need to have your equipment inspected and checked for faults or potential future problems. Just get in touch to discuss your equipment and make a booking.

Breakdown Response

Fully prepared, with a team of experienced engineers and an availability of parts for all the couches and chairs in the range, Knight Imaging is ideally placed to diagnose and repair in times of emergency.

Knight Imaging Limited

Paradigm House, 3 Melbourne Avenue
Broomhill, Sheffield, S10 2QJ, UK

E-mail: sales@knightimaging.co.uk

Tel: +44 (0)114 267 0482

Fax: +44 (0)114 266 3944

www.knightimaging.co.uk

